

Anti-TRIM28 Antibody

Catalog Number: A00409


Overview

Product Name	Anti-TRIM28 Antibody
Reactive Species	Human, Mouse, Rat
Description	Boster Bio Anti-TRIM28 Antibody catalog # A00409. Tested in WB,IHC,ICC/IF,IP applications. This antibody reacts with Human,Mouse,Rat.
Application	IP, IF, IHC, ICC, WB
Clonality	Polyclonal
Formulation	Rabbit IgG, 1mg/ml in PBS with 0.02% sodium azide, 50% glycerol, pH7.2
Storage Instructions	Store at -20°C for one year. For short term storage and frequent use, store at 4°C for up to one month. Avoid repeated freeze-thaw cycles.
Host	Rabbit
Uniprot ID	Q13263


Technical Details

Immunogen	Recombinant fusion protein of human KAP1/KAP1/TRIM28(NP_005753.1).
Cross Reactivity	Does not cross-react with eNOS or nNOS.
Isotype	IgG
Form	Liquid
Concentration	1 mg/ml
Purification	The antibody was affinity-purified from rabbit antiserum by affinity-chromatography using epitope-specific immunogen and the purity is > 95% (by SDS-PAGE).
Suggested Dilutions	WB: 1:1000-1:3000 IHC: 1:100-1:200 ICC/IF: 1:50-1:200 IP: 1:50-1:200


Anti-TRIM28 Antibody (A00409) Images

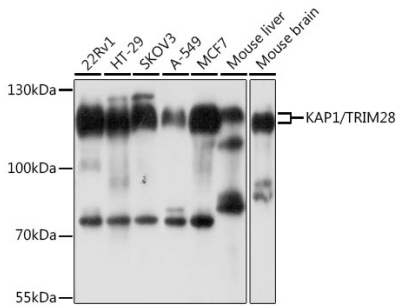


Western blot analysis of extracts of various cell lines, using KAP1/KAP1/TRIM28 antibody at 1:1000 dilution. Secondary antibody: HRP Goat Anti-Rabbit IgG at 1:10000 dilution. Lysates/proteins: 25ug per lane. Blocking buffer: 3% nonfat dry milk in TBST. Detection: ECL Basic Kit. Exposure time: 5s.

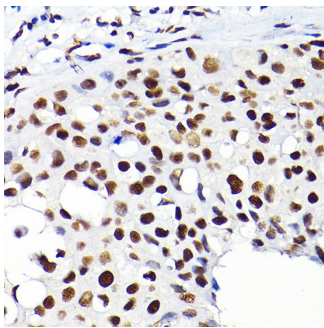

Immunohistochemistry of paraffin-embedded human breast cancer using KAP1/KAP1/TRIM28 Rabbit pAb at dilution of 1:200. Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

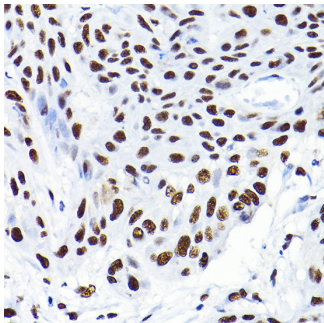

Immunohistochemistry of paraffin-embedded human esophageal cancer using KAP1/KAP1/TRIM28 Rabbit pAb at dilution of 1:200. Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

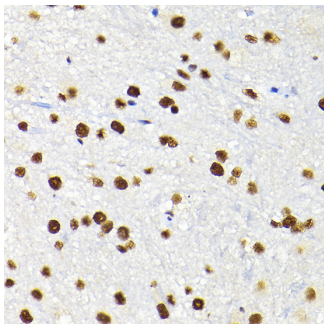

Immunohistochemistry of paraffin-embedded mouse spinal cord using KAP1/KAP1/TRIM28 Rabbit pAb at dilution of 1:200. Perform high pressure antigen retrieval with 10 mM citrate buffer pH 6.0 before commencing with IHC staining protocol.

Submit a product review to [Biocompare.com](https://www.biocompare.com)

Submit a review of this product to [Biocompare.com](https://www.biocompare.com) to receive a \$20 Amazon.com giftcard! Your reviews help your fellow scientists make the right decisions. Thank you for your contribution.



Anti-TRIM28 Antibody

For Research Use Only. Not for use in diagnostic procedures.