

888-466-3604 | support@bosterbio.com | www.bosterbio.com

Anti-AMPK alpha 1/PRKAA1 Antibody Picoband™

Catalog Number: A00994-3

About PRKAA1

5'-AMP-activated protein kinase catalytic subunit alpha-1 is an enzyme that in humans is encoded by the PRKAA1 gene. The protein encoded by this gene belongs to the ser/thr protein kinase family. It is the catalytic subunit of the 5'-prime-AMP-activated protein kinase (AMPK). AMPK is a cellular energy sensor conserved in all eukaryotic cells. The kinase activity of AMPK is activated by the stimuli that increase the cellular AMP/ATP ratio. AMPK regulates the activities of a number of key metabolic enzymes through phosphorylation. It protects cells from stresses that cause ATP depletion by switching off ATP-consuming biosynthetic pathways.

Overview

Product Name	Anti-AMPK alpha 1/PRKAA1 Antibody Picoband™
Reactive Species	Human, Mouse, Rat
Description	Boster Bio Anti-AMPK alpha 1/PRKAA1 Antibody Picoband [™] catalog # A00994-3. Tested in ELISA, Flow Cytometry, IF, ICC, WB applications. This antibody reacts with Human, Mouse, Rat.
Application	ELISA, Flow Cytometry, IF, ICC, WB
Clonality	Polyclonal
Formulation	Each vial contains 4mg Trehalose, 0.9mg NaCl, 0.2mg Na $_2$ HPO $_4$, 0.05mg NaN $_3$.
Storage Instructions	Store at -20°C for one year from date of receipt. After reconstitution, at 4°C for one month. It can also be aliquotted and stored frozen at -20°C for six months. Avoid repeated freeze-thaw cycles.
Host	Rabbit
Uniprot ID	Q13131

Technical Details

Immunogen	E. coli-derived human AMPK alpha 1 recombinant protein (Position: F377-R446).
Predicted Reactive Species	Human
Recommended Detection Systems	Boster recommends Enhanced Chemiluminescent Kit with anti-Rabbit IgG (EK1002) for Western blot, and HRP Conjugated anti-Rabbit IgG Super Vision Assay Kit (SV0002-1) for ICC.
Cross Reactivity	No cross-reactivity with other proteins.
Isotype	Rabbit IgG
Form	Lyophilized
Concentration	Adding 0.2 ml of distilled water will yield a concentration of 500 ug/ml.



BOSTER BIOLOGICAL TECHNOLOGY 3942 B Valley Ave, Pleasanton, CA 94566

888-466-3604 | support@bosterbio.com | www.bosterbio.com

Purification	Immunogen affinity purified.
Suggested Dilutions	Dilute the sample so that the expected range of concentrations fall within the detection range of this kit. If the expected range of concentration is unknown, a pilot test should be conducted to decide the optimal dilution ratio for your samples. Some PubMed article(s) citing the expression level of this target are as follows: Boster Bio's internal QC testing used: Western blot, 0.1-0.5ug/ml Immunocytochemistry/Immunofluorescence, 2ug/ml Flow Cytometry, 1-3ug/1x10 ⁶ cells Direct ELISA, 0.1-0.5ug/ml



888-466-3604 | support@bosterbio.com | www.bosterbio.com

Anti-AMPK alpha 1/PRKAA1 Antibody Picoband[™] (A00994-3) Images

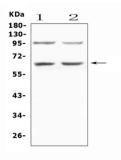


Figure 1. Western blot analysis of AMPK alpha 1 using anti-AMPK alpha 1 antibody (A00994-3).

Electrophoresis was performed on a 5-20% SDS-PAGE gel at 70V (Stacking gel) / 90V (Resolving gel) for 2-3 hours. The sample well of each lane was loaded with 50ug of sample under reducing conditions. Lane 1: human PANC-1 whole cell lysates, Lane 2: human lurkat whole cell lysates.

After Electrophoresis, proteins were transferred to a Nitrocellulose membrane at 150mA for 50-90 minutes. Blocked the membrane with 5% Non-fat Milk/ TBS for 1.5 hour at RT. The membrane was incubated with rabbit anti-AMPK alpha 1 antigen affinity purified polyclonal antibody (Catalog # A00994-3) at 0.5 ug/mL overnight at 4°C, then washed with TBS-0.1%Tween 3 times with 5 minutes each and probed with a goat anti-rabbit IgG-HRP secondary antibody at a dilution of 1:10000 for 1.5 hour at RT. The signal is developed using an Enhanced Chemiluminescent detection (ECL) kit (Catalog # EK1002) with Tanon 5200 system. A specific band was detected for AMPK alpha 1 at approximately 63KD. The expected band size for AMPK alpha 1 is at 63KD.

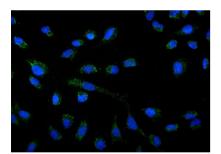


Figure 2. IF analysis of AMPK alpha 1 using anti-AMPK alpha 1 antibody (A00994-3).

AMPK alpha 1 was detected in immunocytochemical section of A549 cells. Enzyme antigen retrieval was performed using IHC enzyme antigen retrieval reagent (AR0022) for 15 mins. The cells were blocked with 10% goat serum. And then incubated with 2ug/mL rabbit anti-AMPK alpha 1 Antibody (A00994-3) overnight at 4°C. DyLight®488 Conjugated Goat Anti-Rabbit IgG (BA1127) was used as secondary antibody at 1:100 dilution and incubated for 30 minutes at 37°C. The section was counterstained with DAPI. Visualize using a fluorescence microscope and filter sets appropriate for the label used.

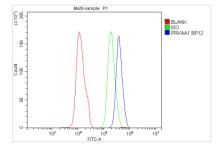


Figure 3. Flow Cytometry analysis of SiHa cells using anti-AMPK alpha 1 antibody (A00994-3).

Overlay histogram showing SiHa cells stained with A00994-3 (Blue line).The cells were blocked with 10% normal goat serum. And then incubated with rabbit anti-AMPK alpha 1 Antibody (A00994-3,1ug/1x10⁶ cells) for 30 min at 20°C. DyLight®488 conjugated goat anti-rabbit IgG (BA1127, 5-10ug/1x10⁶ cells) was used as secondary antibody for 30 minutes at 20°C. Isotype control antibody (Green line) was rabbit IgG (1ug/1x10⁶) used under the same conditions. Unlabelled sample (Red line) was also used as a control.

2 Publications Citing This Product



888-466-3604 | support@bosterbio.com | www.bosterbio.com

Protective Effect Of Glycyrrhizic Acid On Renal Tubular Epithelial Cell Injury Induced By High Glucose.

2. PubMed ID: 29033843, The apolipoprotein A?I mimetic peptide, D?4F, restrains neointimal formation through heme oxygenase?1 up?regulation

Visit bosterbio.com/anti-ampk-alpha-1-antibody-a00994-3-boster.html to see all 2 publications.

Submit a product review to Biocompare.com



Submit a review of this product to Biocompare.com to receive a \$20 Amazon.com giftcard! Your reviews help your fellow scientists make the right decisions. Thank you for your contribution.

Anti-AMPK alpha 1/PRKAA1 Antibody ™