

Anti-Psoriasin/S100A7 Antibody Picoband™

Catalog Number: A02369-1

About S100A7

S100 calcium-binding protein A7 (S100A7), also known as psoriasin, is a protein that in humans is encoded by the S100A7 gene. The protein encoded by this gene is a member of the S100 family of proteins containing 2 EF-hand calcium-binding motifs. S100 proteins are localized in the cytoplasm and/or nucleus of a wide range of cells, and involved in the regulation of a number of cellular processes such as cell cycle progression and differentiation. S100 genes include at least 13 members which are located as a cluster on chromosome 1q21. This protein differs from the other S100 proteins of known structure in its lack of calcium binding ability in one EF-hand at the N-terminus. The protein is overexpressed in hyperproliferative skin diseases, exhibits antimicrobial activities against bacteria and induces immunomodulatory activities.

Overview

| Product Name | Anti-Psoriasin/S100A7 Antibody Picoband™ |
|----------------------|---|
| Reactive Species | Human, Mouse, Rat |
| Description | Boster Bio Anti-Psoriasin/S100A7 Antibody Picoband™ catalog # A02369-1. Tested in ELISA, Flow Cytometry, IF, IHC, ICC, WB applications. This antibody reacts with Human, Mouse, Rat. |
| Application | ELISA, Flow Cytometry, IF, IHC, ICC, WB |
| Clonality | Polyclonal |
| Formulation | Each vial contains 4mg Trehalose, 0.9mg NaCl, 0.2mg Na ₂ HPO ₄ , 0.05mg NaN ₃ . |
| Storage Instructions | Store at -20°C for one year from date of receipt. After reconstitution, at 4°C for one month. It can also be aliquotted and stored frozen at -20°C for six months. Avoid repeated freeze-thaw cycles. |
| Host | Rabbit |
| Uniprot ID | P31151 |

Technical Details

| Immunogen | E. coli-derived human Psoriasin recombinant protein (Position: S2-Q101). |
|-------------------------------|--|
| Recommended Detection Systems | Boster recommends Enhanced Chemiluminescent Kit with anti-Rabbit IgG (EK1002) for Western blot, and HRP Conjugated anti-Rabbit IgG Super Vision Assay Kit (SV0002-1) for IHC(P) and ICC. |
| Cross Reactivity | No cross-reactivity with other proteins. |
| Isotype | Rabbit IgG |
| Form | Lyophilized |
| Concentration | Adding 0.2 ml of distilled water will yield a concentration of 500 ug/ml. |
| | |





888-466-3604 | support@bosterbio.com | www.bosterbio.com

| Purification | Immunogen affinity purified. |
|---------------------|---|
| Suggested Dilutions | Dilute the sample so that the expected range of concentrations fall within the detection range of this kit. If the expected range of concentration is unknown, a pilot test should be conducted to decide the optimal dilution ratio for your samples. Some PubMed article(s) citing the expression level of this target are as follows: Boster Bio's internal QC testing used: Western blot, 0.1-0.5ug/ml Immunohistochemistry (Paraffin-embedded Section), 0.5-1ug/ml Immunocytochemistry/Immunofluorescence, 2ug/ml Flow Cytometry, 1-3ug/1x10 ⁶ cells Direct ELISA, 0.1-0.5ug/ml |



Anti-Psoriasin/S100A7 Antibody Picoband™ (A02369-1) Images

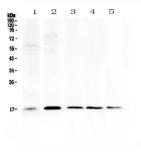


Figure 1. Western blot analysis of Psoriasin using anti-Psoriasin antibody (A02369-1).

Electrophoresis was performed on a 5-20% SDS-PAGE gel at 70V (Stacking gel) / 90V (Resolving gel) for 2-3 hours. The sample well of each lane was loaded with 50ug of sample under reducing conditions.

Lane 1: human placenta tissue lysates,

Lane 2: human Hela cell lysates,

Lane 3: rat spleen tissue lysates,

Lane 4: rat RH35 cell lysates,

Lane 5: mouse spleen tissue lysates.

After Electrophoresis, proteins were transferred to a Nitrocellulose membrane at 150mA for 50-90 minutes. Blocked the membrane with 5% Non-fat Milk/ TBS for 1.5 hour at RT. The membrane was incubated with rabbit anti-Psoriasin antigen affinity purified polyclonal antibody (Catalog # A02369-1) at 0.5 ug/mL overnight at 4°C, then washed with TBS-0.1%Tween 3 times with 5 minutes each and probed with a goat anti-rabbit IgG-HRP secondary antibody at a dilution of 1:10000 for 1.5 hour at RT. The signal is developed using an Enhanced Chemiluminescent detection (ECL) kit (Catalog # EK1002) with Tanon 5200 system. A specific band was detected for Psoriasin at approximately 17KD. The expected band size for Psoriasin is at 11KD.

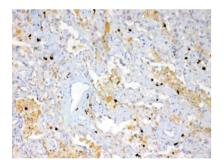


Figure 2. IHC analysis of Psoriasin using anti-Psoriasin antibody (A02369-1).

Psoriasin was detected in paraffin-embedded section of human lung cancer tissue. Heat mediated antigen retrieval was performed in citrate buffer (pH6, epitope retrieval solution) for 20 mins. The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 1ug/ml rabbit anti-Psoriasin Antibody (A02369-1) overnight at 4°C. Biotinylated goat anti-rabbit IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using Strepavidin-Biotin-Complex (SABC)(Catalog # SA1022) with DAB as the chromogen.

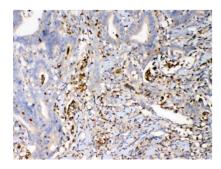


Figure 3. IHC analysis of Psoriasin using anti-Psoriasin antibody (A02369-1).

Psoriasin was detected in paraffin-embedded section of human rectal cancer tissue. Heat mediated antigen retrieval was performed in citrate buffer (pH6, epitope retrieval solution) for 20 mins. The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 1ug/ml rabbit anti-Psoriasin Antibody (A02369-1) overnight at 4°C. Biotinylated goat anti-rabbit IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using Strepavidin-Biotin-Complex (SABC)(Catalog # SA1022) with DAB as the chromogen.



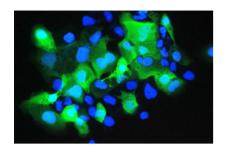


Figure 4. IF analysis of Psoriasin using anti-Psoriasin antibody (A02369-1).

Psoriasin was detected in immunocytochemical section of A431 cells. Enzyme antigen retrieval was performed using IHC enzyme antigen retrieval reagent (AR0022) for 15 mins. The cells were blocked with 10% goat serum. And then incubated with 2ug/mL rabbit anti-Psoriasin Antibody (A02369-1) overnight at 4°C. Cy3 Conjugated Goat Anti-Rabbit IgG (BA1032) was used as secondary antibody at 1:100 dilution and incubated for 30 minutes at 37°C. The section was counterstained with DAPI. Visualize using a fluorescence microscope and filter sets appropriate for the label used.

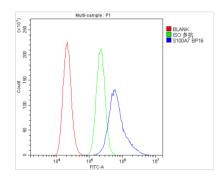


Figure 5. Flow Cytometry analysis of A431 cells using anti-Psoriasin antibody (A02369-1).

Overlay histogram showing A431 cells stained with A02369-1 (Blue line). The cells were blocked with 10% normal goat serum. And then incubated with rabbit anti-Psoriasin Antibody (A02369-1,1ug/1x10 6 cells) for 30 min at 20 $^\circ$ C. DyLight® 488 conjugated goat anti-rabbit IgG (BA1127, 5-10ug/1x10 6 cells) was used as secondary antibody for 30 minutes at 20 $^\circ$ C. Isotype control antibody (Green line) was rabbit IgG (1ug/1x10 6) used under the same conditions. Unlabelled sample (Red line) was also used as a control.

Submit a product review to Biocompare.com





Submit a review of this product to Biocompare.com to receive a \$20 Amazon.com giftcard! Your reviews help your fellow scientists make the right decisions. Thank you for your contribution.

Anti-Psoriasin/S100A7 Antibody ™