

Anti-TDRD5 Antibody Picoband®

Catalog Number: A13514

About TDRD5

Predicted to be involved in P granule organization; spermatid development; and transposable element silencing by piRNA-mediated DNA methylation. Predicted to be located in chromatoid body and pi-body. Predicted to be active in synapse.

Submit a product review to Biocompare.com

[Submit a review of this product](#) to Biocompare.com to receive a \$20 Amazon.com giftcard! Your reviews help your fellow scientists make the right decisions. Thank you for your contribution.



Overview


Product Name	Anti-TDRD5 Antibody Picoband®
Reactive Species	Human, Mouse, Rat
Description	Boster Bio Anti-TDRD5 Antibody Picoband® catalog # A13514. Tested in WB, Flow Cytometry, ELISA applications. This antibody reacts with Human, Mouse, Rat. The brand Picoband indicates this is a premium antibody that guarantees superior quality, high affinity, and strong signals with minimal background in Western blot applications. Only our best-performing antibodies are designated as Picoband, ensuring unmatched performance.
Application	ELISA, Flow Cytometry, WB
Clonality	Polyclonal
Formulation	Each vial contains 4 mg Trehalose, 0.9 mg NaCl, 0.2 mg Na ₂ HPO ₄ .
Storage Instructions	At -20°C for one year from date of receipt. After reconstitution, at 4°C for one month. It can also be aliquotted and stored frozen at -20°C for six months. Avoid repeated freezing and thawing.
Host	Rabbit
Uniprot ID	Q8NAT2

Technical Details


Immunogen	E.coli-derived human TDRD5 recombinant protein (Position: M1-E722). Human TDRD5 shares 82.5% and 82.7% amino acid (aa) sequence identity with mouse and rat TDRD5, respectively.
Form	Lyophilized
Concentration	Adding 0.2 ml of distilled water will yield a concentration of 500 ug/ml.

Purification	Immunogen affinity purified.
Suggested Dilutions	Western blot, 0.25-0.5 ug/ml, Human, Mouse, Rat Flow Cytometry (Fixed), 1-3 ug/1x10 ⁶ cells, Human ELISA, 0.1-0.5 ug/ml Note: For flow cytometry (FCM), validated for fixed cells only. Not validated for FACS or surface flow cytometry using live cells.

Anti-TDRD5 Antibody Picoband® (A13514) Images

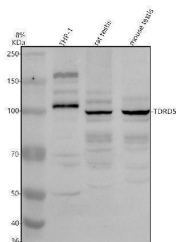

Western blot analysis of TDRD5 using anti-TDRD5 antibody (A13514). Electrophoresis was performed on a 8% SDS-PAGE gel at 80V (Stacking gel) / 120V (Resolving gel) for 2 hours. The sample well of each lane was loaded with 30 ug of sample under reducing conditions. Lane 1: human THP-1 whole cell lysates, Lane 2: rat testis tissue lysates, Lane 3: mouse testis tissue lysates. After electrophoresis, proteins were transferred to a nitrocellulose membrane at 150 mA for 50-90 minutes. Blocked the membrane with 5% non-fat milk/TBS for 1.5 hour at RT. The membrane was incubated with rabbit anti-TDRD5 antigen affinity purified polyclonal antibody (A13514) at 0.5 ug/mL overnight at 4°C, then washed with TBS-0.1%Tween 3 times with 5 minutes each and probed with a goat anti-rabbit IgG-HRP secondary antibody at a dilution of 1:5000 for 1.5 hour at RT. The signal is developed using an ECL Plus Western Blotting Substrate (Catalog # AR1196-200) with Tanon 5200 system. A specific band was detected for TDRD5 at approximately 110-120 kDa. The expected band size for TDRD5 is at 110 kDa.

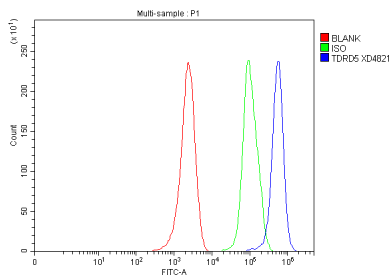

Flow Cytometry analysis of MCF-7 cells using anti-TDRD5 antibody (A13514). Overlay histogram showing MCF-7 cells stained with A13514 (Blue line). To facilitate intracellular staining, cells were fixed with 4% paraformaldehyde and permeabilized with permeabilization buffer. The cells were blocked with 10% normal goat serum. And then incubated with rabbit anti-TDRD5 Antibody (A13514, 1 ug/1x10⁶ cells) for 30 min at 20°C. DyLight® 488 conjugated goat anti-rabbit IgG (BA1127, 5-10 ug/1x10⁶ cells) was used as secondary antibody for 30 minutes at 20°C. Isotype control antibody (Green line) was rabbit IgG (1 ug/1x10⁶) used under the same conditions. Unlabelled sample without incubation with primary antibody and secondary antibody (Red line) was used as a blank control.

Anti-TDRD5 Antibody

For Research Use Only. Not for use in diagnostic procedures.