

PicoKine[™] ELISA

Catalog number: EK0310-BV

For the quantitation of **Bovine BMP5** concentrations in cell culture supernatants, serum and bone tissue.

This package insert must be read in its entirety before using this product. For research use only. Not for use in diagnostic procedures.



Bovine BMP-5 ELISA Kit PicoKine®

Catalog Number: EK0310-BV

Assay Principle

The Boster Picokine[™] Bovine BMP5 Pre-Coated ELISA (Enzyme-Linked Immunosorbent Assay) kit is a solid-phase immunoassay specially designed to measure Bovine BMP5 with a 96-well strip plate that is pre-coated with antibody specific for BMP5. The detection antibody is a biotinylated antibody specific for BMP5. The capture antibody is polyclonal antibody from goat and the detection antibody is polyclonal antibody from goat. The kit includes Bovine BMP5 protein as standards.

To measure Bovine BMP5, add standards and samples to the wells, then add the biotinylated detection antibody. Wash the wells with PBS or TBS buffer, and add Avidin-Biotin-Peroxidase Complex (ABC-HRP). Wash away the unbounded ABC-HRP with PBS or TBS buffer and add TMB. TMB is an HRP substrate and will be catalyzed to produce a blue color product, which changes into yellow after adding the acidic stop solution. The absorbance of the yellow product at 450nm is linearly proportional to Bovine BMP5 in the sample. Read the absorbance of the yellow product in each well using a plate reader, and benchmark the sample wells' readings against the standard curve to determine the concentration of Bovine BMP5 in the sample. For more information on assay principle, protocols, and troubleshooting tips, see Boster's ELISA Resource Center at https://www.bosterbio.com/elisa-technical-resource-center.

Overview

Product Name	Bovine BMP-5 ELISA Kit PicoKine®
Reactive Species	Bovine
Size	96 wells/kit, with removable strips.
Description	Bovine BMP-5 ELISA Kit PicoKine™ (96 Tests). Quantitate Bovine BMP5 in cell culture supernatants, serum and bone tissue. Sensitivity: 10pg/ml.
Sensitivity*	<10 pg/ml
Detection Range	156 pg/ml - 10,000 pg/ml
Storage Instructions	Store at 4°C for 6 months, at -20°C for 12 months. Avoid multiple freeze-thaw cycles (Ships with gel ice, can store for up to 3 days in room temperature. Freeze upon receiving.)
Uniprot ID	P22003

*The sensitivity or the minimum detectable dose (MDD) is the lower limit of the target protein that can be detected by the kit. It is determined by adding two standard deviations to the mean O.D. value of twenty (20) blank wells and calculating the corresponding concentration.

Technical Details

Capture/Detection Antibodies	The capture antibody is polyclonal antibody from goat and the detection antibody is polyclonal antibody from goat.

BOSTER BIOLOGICAL TECHNOLOGY 3942 B Valley Ave, Pleasanton, CA 94566



888-466-3604 | support@bosterbio.com | www.bosterbio.com

Specificity	Natural and recombinant Bovine BMP5
Standard	Expression system for standard: NS0; Immunogen sequence: A317-H454
Cross-reactivity	This kit is for the detection of Bovine BMP5. No significant cross-reactivity or interference between BMP5 and its analogs was observed. This claim is limited by existing techniques; therefore, cross-reactivity may exist with untested analogs.

Preparations Before Assay

Please read the following instructions before starting the experiment.

1. Read this manual in its entirety in order to minimize the chance of error.

2. Confirm that you have the appropriate non-supplied equipment available.

3. Confirm that the species, target antigen, and sensitivity of this kit are appropriate for your intended application.

4. Confirm that your samples have been prepared appropriately based upon recommendations (see Sample Preparation) and that you have sufficient sample volume for use in the assay.

5. When first using a kit, appropriate validation steps should be taken before using valuable samples. Confirm that the kit adequately detects the target antigen in your intended sample type(s) by running control samples.

6. If the concentration of target antigen within your samples is unknown, a preliminary experiment should be run using a control sample to determine the optimal sample dilution (see Sample Preparation).

7. To inspect the validity of experiment operation and the appropriateness of sample dilution proportion, a pilot experiment using standards and a small number of samples is recommended.

8. Before using the kit, spin tubes to bring down all components to the bottom of the tubes.

9. Don't let the 96-well plate dry out since this will inactivate active components on the plate.

10. Don't reuse tips and tubes to avoid cross-contamination.

11. Avoid using the reagents from different batches together.

12. The kit should not be used beyond the expiration date on the kit label. Any variation in diluent, operator, pipetting technique, washing technique, incubation time or temperature, and kit age can cause variation in binding. Variations in sample collection, processing, and storage may cause sample value differences.



Kit Components/Materials Provided

Description	Quantity	Volume	Storage of opened/reconstituted material	
Anti-Bovine BMP5 Pre-coated 96-well Strip Microplate	1	12 strips of 8 wells	Return unused wells to the foil pouch. Reseal along the entire edge of the zip-seal. May be stored for up to 1 month at 4°C provided this is within the expiration date of the kit.	
Bovine BMP5 Standard	2	10 ng/tube	Discard the BMP5 stock solution after 12 hours at 4°C. May be stored at -20°C for 48 hours.	
Bovine BMP5 Biotinylated Antibody (100x)	1	100 µl	May be stored for up to 1 month at	
Avidin-Biotin-Peroxidase Complex (100x)	1	100 µl	 4°C provided this is within the expiration date of the kit. 	
Sample Diluent	1	30 ml		
Antibody Diluent	1	12 ml		
Avidin-Biotin-Peroxidase Diluent	1	12 ml		
Color Developing Reagent (TMB)	1	10 ml		
Stop Solution	1	10 ml		
Wash Buffer (25x)	1	20 ml		
Plate Sealers	4	Piece		

Required Materials That Are Not Supplied

 $\label{eq:main_star} \text{Microplate reader capable of reading absorbance at 450 \, nm.}$

Incubator.

Automated plate washer (optional)

Pipettes and pipette tips capable of precisely dispensing 0.5 µl through 1 ml volumes of aqueous solutions.

Multichannel pipettes are recommended for a large numbers of samples.

Deionized or distilled water.

500 ml graduated cylinders.

Test tubes for dilution.




Bovine BMP5 ELISA Standard Curve Example

The highest O.D. value might be higher or lower than in the example. The experiment result is statistically significant if the highest O.D. value is no less than 1.0.

Concentratio	on0	156	312	625	125000	25000	5000	10000
(pg/ml) O.D.	0.056	0.099	0.116	0.201	0.348	0.634	1.135	1.946

Bovine BMP-5 PicoKine ELISA Kit standard curve

A standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

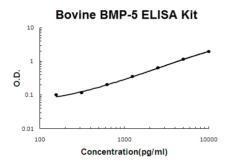

Intra/Inter-Assay Variability

Boster spends great efforts in documenting lot-to-lot variability and ensuring our assay kits produce robust data that are reproducible.

Intra-Assay Precision (Precision within an assay): Three samples of known concentration were tested on one plate to assess intra-assay precision.

Inter-Assay Precision (Precision across assays): Three samples of known concentration were tested in separate assays to assess inter-assay precision.

	Intra-Assay Precision			Inter-Assay Precision		
Sample	1	2	3	1	2	3
n	16	16	16	24	24	24
Mean (ng/ml)	1.23	3.46	6.92	1.38	3.57	6.85
Standard deviation	0.059	0.180	0.388	0.073	0.203	0.425
CV (%)	4.8	5.2	5.6	5.3	5.7	6.2



Reproducibility

We ensure reproducibility by testing three samples with differing concentrations of BMP5 in ELISA kits from four different production batches/lots.

Preparation Before The Experiment

Item	Preparation
All reagents	Bring all reagents to room temperature (18-25°C) prior to use. Please DO NOT equilibrate unused plate well strips to room temperature. They should be sealed and stored in the original packaging. The assay can also be done at room temperature however we recommend doing it at 37°C for best consistency with our QC results. Also, the TMB incubation time estimate (15-25 min) is based on incubation at 37°C.
Wash buffer	Prepare 500 ml of Working Wash Buffer by diluting the supplied 20 ml of Wash Buffer (25 x) with 480 ml of deionized or distilled water. If crystals have formed in the concentrate, warm to room temperature and mix it gently until crystals have completely dissolved.
Biotinylated Anti-Bovine BMP5 antibody	It is recommended to prepare this reagent immediately prior to use by diluting the Bovine BMP5 Biotinylated antibody (100x) 1:100 with Antibody Diluent. Prepare 100 µl by adding 1 µl of Biotinylated antibody (100x) to 99 µl of Antibody Diluent for each well. Mix gently and thoroughly and use within 2 hours of generation.
Avidin-Biotin-Peroxidase Complex	It is recommended to prepare this reagent immediately prior to use by diluting the Avidin-Biotin-Peroxidase Complex (100x) 1:100 with Avidin-Biotin-Peroxidase Diluent. Prepare 100 µl by adding 1 µl of Avidin-Biotin- Peroxidase Complex (100x) to 99 µl of Avidin-Biotin-Peroxidase Diluent for each well. Mix gently and thoroughly and use within 2 hours of generation.
Bovine BMP5 Standard	It is recommended that the standards be prepared no more than 2 hours prior to performing the experiment. Use one 10 ng of lyophilized Bovine BMP5 standard for each experiment. Gently spin the vial prior to use. Reconstitute the standard to a stock concentration of 10 ng/ml using 1ml of sample diluent. Allow the standard to sit for a minimum of 10 minutes with gentle agitation prior to making dilutions.





Microplate	The included microplate is coated with capture antibodies and is ready-to-use. It does not require additional washing or blocking. The unused well strips should be sealed and stored in the original packaging.
Samples	Dilute the sample so that the expected range of concentrations fall within the detection range of this kit. If the expected range of concentration is unknown, a pilot test should be conducted to decide the optimal dilution ratio for your samples. Some PubMed article(s) citing the expression level of this target are as follows:

Dilution of Bovine BMP5 Standard

1. Number tubes 1-8. Final Concentrations to be Tube # 1: 10,000.00 pg/ml, # 2: 5,000.00 pg/ml, # 3: 2,500.00 pg/ml, # 4: 1,250.00 pg/ml,

5: 625.00 pg/ml, # 6: 312.50 pg/ml, # 7: 156.25 pg/ml, # 8: Sample Diluent serves as the zero standard (0 pg/ml).

- 2. For standard #1, add 1000 µl of undiluted standard stock solution to tube #1.
- 3. Add 300 μ l of sample diluent to tubes # 2-7.
- 4. To generate standard # 2, add 300 μl of standard # 1 from tube # 1 to tube # 2 for a final volume of 600 μl. Mix thoroughly.
- 5. To generate standard # 3, add 300 µl of standard # 2 from tube # 2 to tube # 3 for a final volume of 600 µl. Mix thoroughly.
- 6. Continue the serial dilution for tube # 4-7.

Sample Preparation and Storage

These sample collection instructions and storage conditions are intended as a general guideline, and the sample stability has not been evaluated.

Sample dilution ratios should be determined by a pilot study (run a serial dilution of samples and see which dilution ratio results in the idea O.D., near the middle of the standard range). In general, high concentration samples can be dilutioned by 1:100, mid concentration samples 1:10, low concentration samples 1:2 or neat.

Sample Type	Procedure
Cell culture supernatants	Clear sample of particulates by centrifugation, assay immediately, or store samples at -20°C.
Serum	Use a serum separator tube (SST) and allow serum to clot at room temperature for about four hours. Then, centrifuge for 15 min at approximately 1,000 x g. assay immediately or store samples at -20°C.
Bone tissue	Extract demineralized bone samples in 4M Guanidine-HCl and protease inhibitors. Dissolve the final sample in 2M Guanidine-HCl.



Sample Collection Notes

1. Boster recommends that samples are used immediately upon preparation.

2. Avoid repeated freeze/thaw cycles for all samples.

3. In the event that a sample type not listed above is intended to be used with the kit, it is recommended that the customer conduct validation experiments in order to be confident in the results.

4. Due to chemical interference, the use of tissue or cell extraction samples prepared by chemical lysis buffers may result in inaccurate results.

- 5. Due to factors including cell viability, cell number, or sampling time, samples from cell culture supernatant may not be detected by the kit.
- 6. Samples should be brought to room temperature (18-25°C) before performing the assay without the use of extra heating.

7. Sample concentrations should be predicted before being used in the assay. If the sample concentration is not within the range of the standard curve, users must determine the optimal sample dilutions for their particular experiments.

8. Boster is responsible for the quality and performance of the kit components but is NOT responsible for the performance of customer supplied samples used with the kit.

Sample Dilution Guideline

The user needs to estimate the concentration of the target protein in the sample and use an appropriate dilution factor so that the diluted target protein concentration falls in the range of O.D. values of the standard curve. Dilute the sample using provided diluent buffer. Pilot tests using a dilution series of each sample type are necessary. The sample must be mixed thoroughly with Sample Diluent.

Assay Protocol

It is recommended that all reagents and materials be equilibrated to room temperature (18-25°C) prior to the experiment (see Preparation Before The Experiment, if you have missed this information).

- 1. Prepare all reagents and working standards as directed previously.
- 2. Remove excess microplate strips from the plate frame and seal and store them in the original packaging.

3. Add 100 µl of the standard, samples, or control per well. Add 100 µl of the **Sample Diluent** into the zero well. At least two replicates of each standard, sample, or control is recommended.

4. Cover with the plate sealer provided and incubate for 120 minutes at room temperature (or 90 min. at 37 °C).

5. Remove the cover and discard the liquid in the wells into an appropriate waste receptacle. Invert the plate on the benchtop onto a paper towel and tap the plate to gently blot any remaining liquid. It is recommended that the wells are not allowed to completely dry at any time.

6. Add 100 µl of the prepared 1x Biotinylated Anti-Bovine BMP5 antibody to each well.

7. Cover with a plate sealer and incubate for 90 minutes at room temperature (or 60 minutes at 37°C).

8. Wash the plate 3 times with the **1x wash buffer**:

a. Discard the liquid in the wells into an appropriate waste receptacle. Then, invert the plate on the benchtop onto a paper towel and tap the plate to gently blot any remaining liquid. It is recommended that the wells are not allowed to completely dry at any time.

b. Add 300 µl of the 1x wash buffer to each assay well. (For cleaner background incubate for 60 seconds between each wash).

c. Repeat steps a-b 2 additional times.

d. Discard the wash buffer in the wells into an appropriate waste receptacle. Then, invert the plate on the benchtop onto a paper towel and tap the plate to gently blot any remaining liquid.

9. Add 100 µl of the prepared **1x Avidin-Biotin-Peroxidase Complex** into each well. Cover with the plate sealer provided and incubate for 40 minutes at RT (or 30 minutes at 37°C).

10. Wash the plate 5 times with the **1x wash buffer**:

a. Discard the liquid in the wells into an appropriate waste receptacle. Then, invert the plate on the benchtop onto a paper towel and tap the plate to gently blot any remaining liquid. It is recommended that the wells are not allowed to completely dry at any time.

b. Add 300 µl of the 1x wash buffer to each assay well. (For cleaner background incubate for 60 seconds between each wash).

antibody and ELISA experts c. Repeat steps a-b 4 additional times.

STER

d. Discard the wash buffer in the wells into an appropriate waste receptacle. Then, invert the plate on the benchtop onto a paper towel and tap the plate to gently blot any remaining liquid.

11. Add 90 µl of **Color Developing Reagent** to each well. Cover with the plate sealer provided and incubate in the dark for 30 minutes at RT (or 15-25 minutes at 37°C). (The optimal incubation time must be empirically determined. A guideline to look for is blue shading the top four standard wells, while the remaining standards remain clear.)

12. Add 100 μl of ${\mbox{Stop}}$ Solution to each well. The color should immediately change to yellow.

13. Within 30 minutes of stopping the reaction, the O.D. absorbance should be read with a microplate reader at 450nm.

Assay Protocol Notes

1. Solutions: To avoid cross-contamination, change pipette tips between additions of each standard, between sample additions, and between reagent additions. Also, use separate reservoirs for each reagent.

2. Applying Solutions: All solutions should be added to the bottom of the ELISA plate well. Avoid touching the inside wall of the well. Avoid foaming when possible.

3. Assay Timing: The interval between adding samples to the first and last wells should be minimized. Delays will increase the incubation time differential between wells, which will significantly affect the experimental accuracy and repeatability. For each step in the procedure, total dispensing time for addition of reagents or samples should not exceed 10 minutes.

4. Incubation: To prevent evaporation and ensure accurate results, proper adhesion of plate sealers during incubation steps is necessary. Do not allow wells to sit uncovered for extended periods of time between incubation steps. Do not let wells dry out at any time during the assay. Strictly observe the recommended incubation times and temperatures.

5. Washing: Proper washing procedure is critical. Insufficient washing will result in poor precision and falsely elevated absorbance readings. Residual liquid in the reaction wells should be patted dry against absorbent paper during the washing process. Do not put absorbent paper directly into the reaction wells.

6. Controlling Substrate Reaction Time: After the addition of the TMB Substrate, periodically monitor the color development. Stop color development before the color becomes too deep by adding Stop Solution. The excessively strong color will result in inaccurate absorbance readings.
7. Reading: The microplate reader should be preheated and programmed prior to use. Prior to taking O.D. readings, remove any residual liquid or

fingerprints from the underside of the plate and confirm that there are no bubbles in the wells.

8. Reaction Time Control: Control reaction time should be strictly followed as outlined.

9. Stop Solution: The Stop Solution contains an acid, therefore proper precautions should be taken during its use, such as protection of the eyes, hands, face, and clothing.

10. To minimize the external influence on the assay performance, operational procedures and lab conditions (such as room temperature, humidity, incubator temperature) should be strictly controlled. It is also strongly suggested that the whole assay is performed by the same operator from the beginning to the end.



Data Analysis

Boster Bio offers an easy-to-use online ELISA data analysis tool. Try it out at <u>https://www.bosterbio.com/biology-research-tools/elisa-data-analysis-online</u>

To analyze using manual methods, follow the process below:

Average the duplicate readings for each standard, sample, and control. Subtract the average zero standard O.D. reading. It is recommended that a standard curve be created using computer software to generate a four-parameter logistic (4-PL) curve-fit. A free program capable of generating a four-parameter logistic (4-PL) curve-fit can be found online at: www.myassays.com/four-parameter-logistic-curve.assay. Alternatively, plot the mean absorbance for each standard against the concentration. The measured concentration in the sample can be interpolated by using linear regression of each average relative O.D. against the standard curve generated using curve fitting software. This will generate an adequate but less precise fit of the data.

For diluted samples, the concentration reading from the standard curve must be multiplied by the dilution factor.

Background on BMP5

Bone morphogenetic protein 5 is a protein that in humans is encoded by the BMP5 gene.1,2,3 The protein encoded by this gene is member of the TGFbeta superfamily. Bone morphogenetic proteins are known for their ability to induce bone and cartilage development. BMP5 may play a role in certain cancers. Like other BMP's BMP5 is inhibited by chordin and noggin. It is expressed in the trabecular meshwork and optic nerve head and may have a role in the development and normal function. It is also expressed in the lung and liver. BMP-5 is mapped on Bovine chromosome 6.4 The standard used in this kit is recombinant BMP-5 with the monomer molecular mass of 18KDa.

Submit a Product Review to Biocompare.com



Submit a review of this product to Biocompare.com to receive a \$20 Amazon.com gift card! Your reviews help your fellow scientists make the right decisions. Thank you for your contribution.

Bovine BMP-5 ELISA Kit ®