

Anti-Toll-interacting protein Tollip Antibody Picoband®

Catalog Number: PA2005

About TOLLIP

TOLLIP (TOLL-Interacting Protein), is an inhibitory adaptor protein that in humans is encoded by the TOLLIP gene. Lo et al. (2009) stated that the TOLLIP gene maps to chromosome 11. The mouse gene maps to chromosome 7. By Western blot analysis of embryonic kidney cells, Burns et al. (2000) confirmed the binding of TOLLIP to IL1RAP, to a complex of IL1RAP-IL1R1, and to IL18R. Burns et al. (2000) proposed that IL1B stimulation induces aggregation of IL1Rs, recruitment of MYD88 followed by TOLLIP-IRAK complexes, and the phosphorylation of IRAK by MYD88. This leads to the dissociation of TOLLIP from IRAK, which can then transmit the IL1-induced signals.

Overview


Product Name	Anti-Toll-interacting protein Tollip Antibody Picoband®
Reactive Species	Human, Mouse, Rat
Description	Boster Bio Anti-Toll-interacting protein Tollip Antibody catalog # PA2005. Tested in Flow Cytometry, IP, ICC, IHC, IF, WB applications. This antibody reacts with Human, Mouse, Rat. The brand Picoband indicates this is a premium antibody that guarantees superior quality, high affinity, and strong signals with minimal background in Western blot applications. Only our best-performing antibodies are designated as Picoband, ensuring unmatched performance.
Application	Flow Cytometry, IP, IF, IHC, ICC, WB
Clonality	Polyclonal
Formulation	Each vial contains 4 mg Trehalose, 0.9 mg NaCl and 0.2 mg Na ₂ HPO ₄ .
Storage Instructions	Store at -20°C for one year from date of receipt. After reconstitution, at 4°C for one month. It can also be aliquotted and stored frozen at -20°C for six months. Avoid repeated freeze-thaw cycles.
Host	Rabbit
Uniprot ID	Q9H0E2

Technical Details


Immunogen	A synthetic peptide corresponding to a sequence at the C-terminus of human Tollip.
Recommended Detection Systems	Boster recommends Enhanced Chemiluminescent Kit with anti-Rabbit IgG (EK1002) for Western blot, and HRP Conjugated anti-Rabbit IgG Super Vision Assay Kit (SV0002-1) for IHC(P) and ICC.
Cross Reactivity	No cross-reactivity with other proteins
Isotype	Rabbit IgG
Form	Lyophilized
Concentration	Adding 0.2 ml of distilled water will yield a concentration of 500 ug/ml.

Purification	Immunogen affinity purified.
Suggested Dilutions	Western blot, 0.1-0.5ug/ml, Human, Mouse, Rat Immunohistochemistry (Paraffin-embedded Section), 2-5ug/ml, Human, Mouse, Rat Immunocytochemistry/Immunofluorescence, 5 ug/ml, Human Immunoprecipitation, 0.5-2 ug/ml, Human Flow Cytometry (Fixed), 1-3ug/1x10 ⁶ cells, Human


Anti-Toll-interacting protein Tollip Antibody Picoband® (PA2005) Images

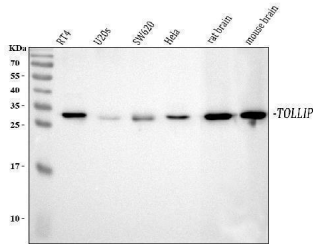

Western blot analysis of Tollip using anti-Tollip antibody (PA2005). Electrophoresis was performed on a 5-20% SDS-PAGE gel at 70V (Stacking gel) / 90V (Resolving gel) for 2-3 hours. The sample well of each lane was loaded with 30 ug of sample under reducing conditions. Lane 1: human RT4 whole cell lysates, Lane 2: human U2OS whole cell lysates, Lane 3: human SW620 whole cell lysates, Lane 4: human Hela whole cell lysates, Lane 5: rat brain tissue lysates, Lane 6: mouse brain tissue lysates. After electrophoresis, proteins were transferred to a nitrocellulose membrane at 150 mA for 50-90 minutes. Blocked the membrane with 5% non-fat milk/TBS for 1.5 hour at RT. The membrane was incubated with rabbit anti-Tollip antigen affinity purified polyclonal antibody (Catalog # PA2005) at 0.5 ug/mL overnight at 4°C, then washed with TBS-0.1%Tween 3 times with 5 minutes each and probed with a goat anti-rabbit IgG-HRP secondary antibody at a dilution of 1:5000 for 1.5 hour at RT. The signal is developed using an Enhanced Chemiluminescent detection (ECL) kit (Catalog # EK1002) with Tanon 5200 system. A specific band was detected for Tollip at approximately 30 kDa. The expected band size for Tollip is at 30 kDa.

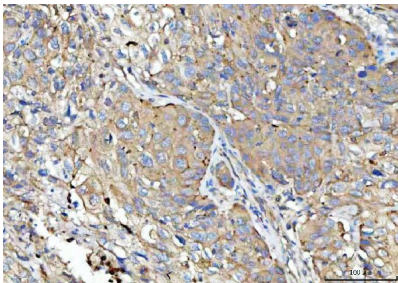


IHC analysis of Tollip using anti-Tollip antibody (PA2005). Tollip was detected in a paraffin-embedded section of human lung cancer tissue. Heat mediated antigen retrieval was performed in EDTA buffer (pH 8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 2 ug/ml rabbit anti-Tollip Antibody (PA2005) overnight at 4°C. Peroxidase Conjugated Goat Anti-rabbit IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using HRP Conjugated Rabbit IgG Super Vision Assay Kit (Catalog # SV0002) with DAB as the chromogen.

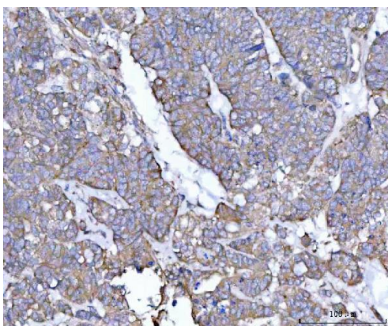


IHC analysis of Tollip using anti-Tollip antibody (PA2005). Tollip was detected in a paraffin-embedded section of human ovarian cancer tissue. Heat mediated antigen retrieval was performed in EDTA buffer (pH 8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 2 ug/ml rabbit anti-Tollip Antibody (PA2005) overnight at 4°C. Peroxidase Conjugated Goat Anti-rabbit IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using HRP Conjugated Rabbit IgG Super Vision Assay Kit (Catalog # SV0002) with DAB as the chromogen.


IHC analysis of Tollip using anti-Tollip antibody (PA2005). Tollip was detected in a paraffin-embedded section of mouse brain tissue. Heat mediated antigen retrieval was performed

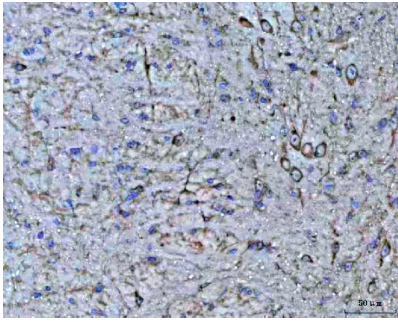


in EDTA buffer (pH 8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 2 ug/ml rabbit anti-Tollip Antibody (PA2005) overnight at 4°C. Peroxidase Conjugated Goat Anti-rabbit IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using HRP Conjugated Rabbit IgG Super Vision Assay Kit (Catalog # SV0002) with DAB as the chromogen.

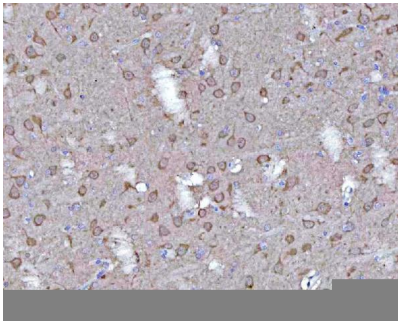

IHC analysis of Tollip using anti-Tollip antibody (PA2005). Tollip was detected in a paraffin-embedded section of rat brain tissue. Heat mediated antigen retrieval was performed in EDTA buffer (pH 8.0, epitope retrieval solution). The tissue section was blocked with 10% goat serum. The tissue section was then incubated with 2 ug/ml rabbit anti-Tollip Antibody (PA2005) overnight at 4°C. Peroxidase Conjugated Goat Anti-rabbit IgG was used as secondary antibody and incubated for 30 minutes at 37°C. The tissue section was developed using HRP Conjugated Rabbit IgG Super Vision Assay Kit (Catalog # SV0002) with DAB as the chromogen.

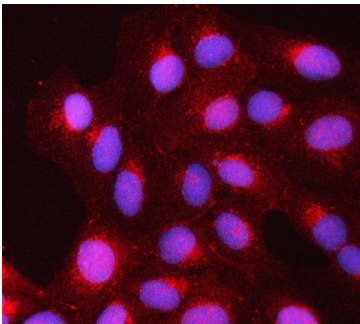


IF analysis of Tollip using anti-Tollip antibody (PA2005). Tollip was detected in an immunocytochemical section of U2OS cells. Enzyme antigen retrieval was performed using IHC enzyme antigen retrieval reagent (AR0022) for 15 mins. The cells were blocked with 10% goat serum. And then incubated with 5 ug/mL rabbit anti-Tollip Antibody (PA2005) overnight at 4°C. Cy3 Conjugated Goat Anti-Rabbit IgG (BA1032) was used as secondary antibody at 1:500 dilution and incubated for 30 minutes at 37°C. The section was counterstained with DAPI. Visualize using a fluorescence microscope and filter sets appropriate for the label used.

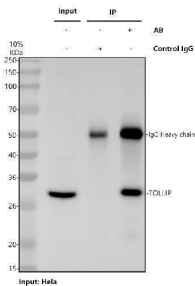

Immunoprecipitating (IP) Tollip in HeLa whole cell lysate. Western blot analysis of Tollip using anti-Tollip antibody (PA2005); Lane 1: HeLa whole cell lysates (30ug); Lane 2: Rabbit control IgG instead of anti-Tollip antibody in HeLa whole cell lysate; Lane 3: anti-Tollip antibody (2ug) + HeLa whole cell lysate (500ug). After electrophoresis, proteins were transferred to a membrane. Then the membrane was incubated with rabbit anti-Tollip antigen affinity purified polyclonal antibody (PA2005) at a dilution of 0.5 ug/mL and probed with a goat anti-rabbit IgG-HRP secondary antibody (Catalog # BA1054). The signal is developed using ECL Plus Western Blotting Substrate (Catalog # AR1196-200). A specific band was detected for Tollip at approximately 30 kDa. The expected band size for Tollip is at 30 kDa.

Flow Cytometry analysis of Daudi cells using anti-Tollip antibody (PA2005). Overlay histogram showing Daudi cells stained with PA2005 (Blue line). To facilitate intracellular staining, cells were fixed with 4% paraformaldehyde and permeabilized with permeabilization buffer. The cells were

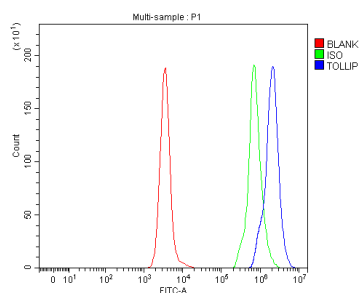

blocked with 10% normal goat serum. And then incubated with rabbit anti-Tollip Antibody (PA2005, 1 ug/1x10⁶ cells) for 30 min at 20°C. DyLight®488 conjugated goat anti-rabbit IgG (BA1127, 5-10 ug/1x10⁶ cells) was used as secondary antibody for 30 minutes at 20°C. Isotype control antibody (Green line) was rabbit IgG (1 ug/1x10⁶) used under the same conditions. Unlabelled sample without incubation with primary antibody and secondary antibody (Red line) was used as a blank control.

1 Publications Citing This Product

1. PubMed ID: 25269519, Liu Mw, Wang Yh, Qian Cy, Li H. Int J Mol Med. 2014 Dec;34(6):1492-504. Doi: 10.3892/Ijmm.2014.1943. Epub 2014 Sep 23. Xuebijing Exerts Protective Effects On Lung Permeability Leakage And Lung Injury By Upregulating Toll-Interacting Protein Expres...

Visit bosterbio.com/anti-tollip-antibody-pa2005-boster.html to see all 1 publications.

Submit a product review to Biocompare.com

Submit a review of this product to Biocompare.com to receive a \$20 Amazon.com giftcard! Your reviews help your fellow scientists make the right decisions. Thank you for your contribution.



Anti-Toll-interacting protein Tollip Antibody

For Research Use Only. Not for use in diagnostic procedures.