



# CRE Luciferase Reporter-HEK293 Cell Line

Catalog number: RC1034

This package insert must be read in its entirety before using this product.  
For research use only. Not for use in diagnostic procedures.

## CRE Luciferase Reporter-HEK293 Cell Line

**Catalog Number:** RC1034, **Storage:** Immediately upon receipt, store in liquid nitrogen. (Ship on dry ice.)

**Contents:** Each vial contains  $2 \sim 3 \times 10^6$  cells in 1 ml of 90% FBS + 10% DMSO.

**Description:** The CRE Luciferase **Reporter cell line** is a stably transfected HEK 293 cell line which expresses Renilla luciferase reporter gene under the transcriptional control of the cAMP response element (CRE). The CRE cell line is designed to monitor the cAMP/PKA signaling pathways and can be used for studying GPCR-linked cAMP/PKA signaling pathways as well as screening of agonists, antagonists or signaling inhibitors related with the cAMP/PKA signaling pathways. Functional activity of the cell line has been validated by serum stimulation assay (Figure 1).

**Applications:** Functional Assay

**Application Notes:** Functional Assay, detecting the transcriptional activity of CRE

**Application Details:**

## Application:

Monitor the GPCR-linked cAMP/PKA signaling pathway. Screen for activators or inhibitors of the GPCR-linked cAMP/PKA signaling pathway.

## Culture conditions:

Cells should be grown at 37°C with 5% CO<sub>2</sub> using DMEM medium supplemented with 10% FBS and 1% Pen/Strep, plus 3 µg/ml of Puromycin. It is recommended to quickly thaw the frozen cells upon receipt or from liquid nitrogen in a 37°C water-bath, transfer to a tube containing 10 ml of growth medium without Puromycin, spin down cells, resuspend cells in pre-warmed growth medium without Puromycin, transfer resuspended cells to T25 flask and culture in 37°C-CO<sub>2</sub> incubator. Leave the T25 flask in the incubator for 2~4 days without disturbing or changing the medium until cells completely recover viability and become adherent. Once cells are over 90% adherent, remove growth medium and passage the cells through trypsinization and centrifugation. At first passage, switch to growth medium containing Puromycin. Cells should be split before they reach complete confluence. To passage the cells, detach cells from culture vessel with Trypsin/EDTA, add complete growth medium and transfer to a tube, spin down cells, resuspend cells and seed appropriate aliquots of cells suspension into new culture vessels. Subcultivation ration = 1:10 to 1:20 weekly.

## Functional validation:

A. Response of CRE HEK293 cells to forskolin. 1. Harvest CRE HEK293 cells and seed cells into a white solid-bottom 96-well microplate in 100 µl of growth medium at  $5 \times 10^4$  c

## CRE Luciferase Reporter-HEK293 Cell Line (RC1034) Images

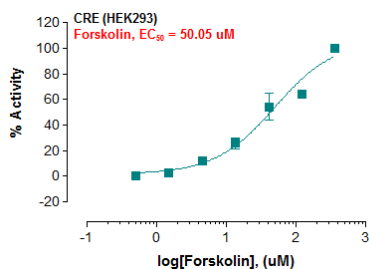


Fig-1: Induction of CRE activity by forskolin in CRE HEK293 cells.

## Submit a product review to Biocompare.com

Submit a review of this product to Biocompare.com to receive a \$20 Amazon giftcard! Your reviews help your fellow scientists make the right decisions. Thank you for your contribution.



CRE Luciferase Reporter-HEK293 Cell Line

For Research Use Only. Not for use in diagnostic procedures.