



TLR5/NF- κ B Luciferase Reporter- HEK293 Cell Line

Catalog number: RC1039

This package insert must be read in its entirety before using this product.
For research use only. Not for use in diagnostic procedures.

TLR5/NF- κ B Luciferase Reporter-HEK293 Cell Line

Catalog Number: RC1039, **Storage:** Immediately upon receipt, store in liquid nitrogen. (Ship on dry ice.)

Contents: Each vial contains $2 \sim 3 \times 10^6$ cells in 1 ml of 90% FBS + 10% DMSO.

Description: The TLR5/NF- κ B Luciferase Reporter cell line is a stably transfected HEK 293 cell line which expresses full-length human Toll-like receptor 5 (TLR5) and Renilla luciferase reporter gene under the transcriptional control of the NF- κ B response element. TLR5 is one of the key innate immune receptors. Functional activity of the cell line has been validated by TLR5 ligand assay, in which upon activation by Flagellin, TLR5 quickly initiates downstream signaling pathway and mediates nuclear translocation of NF- κ B (Figure 1).

Applications: Functional Assay

Application Notes: Functional Assay, detecting the transcriptional activity of TLR5/NF- κ B

Application Details:

Application:

Monitor the TLR5 signaling pathway activity. Screen for activators or inhibitors of the TLR5 signaling pathway.

Culture conditions:

Cells should be grown at 37°C with 5% CO₂ using DMEM medium supplemented with 10% FBS and 1% Pen/Strep, plus 2 μ g/ml of Puromycin and 5 μ g/ml of Blastidicin. It is recommended to quickly thaw the frozen cells upon receipt or from liquid nitrogen in a 37°C water-bath, transfer to a tube containing 10 ml of growth medium without Puromycin and Blastidicin, spin down cells, resuspend cells in pre-warmed growth medium without Puromycin and Blastidicin, transfer resuspended cells to T25 flask and culture in 37°C-CO₂ incubator. Leave the T25 flask in the incubator for 2~4 days without disturbing or changing the medium until cells completely recover viability and become adherent. Once cells are over 90% adherent, remove growth medium and passage the cells through trypsinization and centrifugation. At first passage, switch to growth medium containing Puromycin and Blastidicin. Cells should be split before they reach complete confluence. To passage the cells, detach cells from culture vessel with Trypsin/EDTA, add complete growth medium and transfer to a tube, spin down cells, resuspend cells and seed appropriate aliquots of cells suspension into new culture vessels. Subcultivation ration = 1:10 to 1:20 weekly.

Functional validation:

A. Response of TLR5/NF- κ B - HEK293 cells to Flagellin. 1. Harvest TLR5/NF- κ B - HEK293 cells and seed cells into a white solid-bottom 96-well microplate in 100 μ l of growth medium at 5×10^4 cells/well. 2. Incubate cells at 37°C in a CO₂ incubator for overnight. 3. The next day, stimulate cells with various concentrations of Flagellin. 4. Incubate at 37°C in a CO₂ incubator for 6-16 hours. 5. Add 50 μ l of luciferase assay reagent per well. 6. Incubate at room temperature for 1-5 minutes and measure luminescence using a microplate luminometer.

TLR5/NF- κ B Luciferase Reporter-HEK293 Cell Line (RC1039) Images

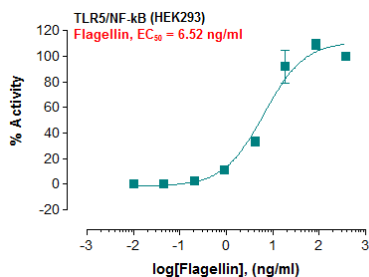


Fig-1: Induction of TLR5 activity by Flagellin in TLR5/NF- κ B HEK293 cells.

Submit a product review to Biocompare.com

Submit a review of this product to Biocompare.com to receive a \$20 Amazon giftcard! Your reviews help your fellow scientists make the right decisions. Thank you for your contribution.



TLR5/NF-kB Luciferase Reporter-HEK293 Cell Line

For Research Use Only. Not for use in diagnostic procedures.