



PicoKine™ ELISA

Catalog number: EK0941

For the quantitation of **Human MIA** concentrations in cell culture supernates, serum, plasma(heparin, EDTA) and saliva.

This package insert must be read in its entirety before using this product. For research use only. Not for use in diagnostic procedures.

Human MIA PicoKine™ ELISA Kit

Catalog Number: EK0941

Assay Principle

The Boster Picokine™ Human MIA Pre-Coated ELISA (Enzyme-Linked Immunosorbent Assay) kit is a solid phase immunoassay specially designed to measure Human MIA with a 96-well strip plate that is pre-coated with antibody specific for MIA. The detection antibody is a biotinylated antibody specific for MIA. The capture antibody is monoclonal antibody from mouse, the detection antibody is polyclonal antibody from goat. The kit contains recombinant Human MIA with immunogen: Expression system for standard: E.coli; Immunogen sequence: G25-131Q. The kit is analytically validated with ready to use reagents.

To measure Human MIA, add standards and samples to the wells, then add the biotinylated detection antibody. Wash the wells with PBS or TBS buffer, and add Avidin-Biotin-Peroxidase Complex (ABC-HRP). Wash away the unbound ABC-HRP with PBS or TBS buffer and add TMB. TMB is substrate to HRP and will be catalyzed to produce a blue color product, which changes into yellow after adding acidic stop solution. The density of the yellow product is linearly proportional to Human MIA in the sample. Read the density of the yellow product in each well using a plate reader, and benchmark the sample wells' readings against the standard curve to determine the concentration of Human MIA in the sample. For more information on assay principle, protocols, and troubleshooting tips, see Boster's ELISA Resource Center at <https://www.bosterbio.com/elisa-technical-resource-center>.

Overview

| | |
|-----------------------------|---|
| Product Name | Human MIA PicoKine™ ELISA Kit |
| Reactive Species | Human |
| Size | 96wells/kit, with removable strips. |
| Description | Sandwich High Sensitivity ELISA kit for Quantitative Detection of Human MIA. 96wells/kit, with removable strips. |
| Sensitivity | <10pg/ml *The sensitivity or the minimum detectable dose (MDD) is the lower limit of target protein that can be detected by the kit. It is determined by adding two standard deviations to the mean O.D. value of twenty (20) blank wells and calculating the corresponding concentration. |
| Detection Range | 15.6pg/ml-1000pg/ml |
| Storage Instructions | Store at 4°C for 6 months, at -20°C for 12 months. Avoid multiple freeze-thaw cycles(Shipped with wet ice.) |
| Uniprot ID | Q16674 |

Technical Details

| | |
|------------------------------|--|
| Capture/Detection Antibodies | The capture antibody is monoclonal antibody from mouse, the detection antibody is polyclonal antibody from goat. |
| Specificity | Natural and recombinant Human MIA |
| Immunogen | Expression system for standard: E.coli; Immunogen sequence: G25-131Q |
| Cross Reactivity | There is no detectable cross-reactivity with other relevant proteins. |

Notice Before Application

Please read the following instructions before starting the experiment.

1. To inspect the validity of experiment operation and the appropriateness of sample dilution proportion, pilot experiment using standards and a small number of samples is recommended.
2. Before using the Kit, spin tubes and bring down all components to the bottom of tubes.
3. Don't let 96-well plate dry, for dry plate will inactivate active components on plate.
4. Don't reuse tips and tubes to avoid cross contamination.
5. Avoid using the reagents from different batches together.

Kit Components/Materials Provided

| Description | Quantity | Volume |
|--|----------|----------------------|
| Anti-Human MIA Pre-coated 96-well strip microplate | 1 | 12 strips of 8 wells |
| Human MIA Standard | 2 | 1ng/tube |
| Human MIA Biotinylated antibody (100x) | 1 | 130 µl |
| Avidin-Biotin-Peroxidase Complex (100x) | 1 | 130 µl |
| Sample Diluent | 1 | 30ml |
| Antibody Diluent | 1 | 12ml |
| Avidin-Biotin-Peroxidase Diluent | 1 | 12ml |
| Color Developing Reagent (TMB) | 1 | 10ml |
| Stop Solution | 1 | 10ml |
| Wash Buffer Powder | 1 | Pack |
| Plate Sealers | 4 | Piece |

Required Materials That Are Not Supplied

Microplate Reader capable of reading absorbance at 450nm.

Automated plate washer (optional)

Pipettes and pipette tips capable of precisely dispensing 0.5 µl through 1 ml volumes of aqueous solutions.

Multichannel pipettes are recommended for large amount of samples.

Deionized or distilled water.

500ml graduated cylinders.


Test tubes for dilution.

Human MIA PicoKine™ ELISA Kit (EK0941) Standard Curve Example

Highest O.D. value might be higher or lower than in the example. The experiment result is statistically significant if the highest O.D. value is no less than 1.0.

| | | | | | | | |
|----------------------------|-------|-------|-------|-------|-------|-------|-------|
| Concentration 0 (pg/ml) | 15.6 | 31.2 | 62.5 | 125 | 250 | 500 | 1000 |
| O.D. | 0.021 | 0.054 | 0.089 | 0.156 | 0.352 | 0.892 | 2.146 |

Human MIA PicoKine ELISA Kit standard curve

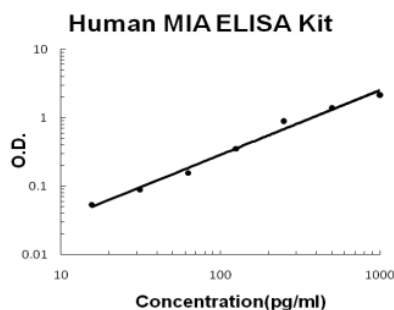

A standard curve is provided for demonstration only. A standard curve should be generated for each set of samples assayed.

Intra/Inter Assay Variability

Boster spend great efforts in documenting lot to lot variability and make sure our assay kits produce robust data that are reproducible.

Intra-Assay Precision (Precision within an assay): Three samples of known concentration were tested on one plate to assess intra-assay precision.

Inter-Assay Precision (Precision across assays): Three samples of known concentration were tested in separate assays to assess inter-assay precision.

| Sample | Intra-Assay Precision | | | Inter-Assay Precision | | |
|--------------------|-----------------------|-------|-------|-----------------------|-------|-------|
| | 1 | 2 | 3 | 1 | 2 | 3 |
| n | 16 | 16 | 16 | 24 | 24 | 24 |
| Mean(pg/ml) | 29 | 172 | 485 | 26 | 165 | 509 |
| Standard deviation | 1.91 | 13.41 | 32.98 | 1.97 | 15.84 | 36.92 |
| CV(%) | 6.6% | 7.8% | 6.8% | 7.6% | 9.6% | 6.9% |

Reproducibility

To assay reproducibility, three samples with differing target protein concentrations were assayed using four different lots.

| Lots | Lot1 (pg/ml) | Lot2 (pg/ml) | Lot3 (pg/ml) | Lot4 (pg/ml) | Mean (pg/ml) | Standard Deviation | CV (%) |
|----------|--------------|--------------|--------------|--------------|--------------|--------------------|--------|
| Sample 1 | 29 | 28 | 27 | 29 | 28 | 0.82 | 2.9% |
| Sample 2 | 172 | 166 | 171 | 170 | 169 | 2.27 | 1.3% |
| Sample 3 | 485 | 453 | 511 | 458 | 476 | 23.22 | 4.8% |

*number of samples for each test n=16.

Preparation Before The Experiment

| Item | Preparation |
|--------------------------------------|--|
| All reagents | Bring all reagents to 37°C prior to use. The assay can also be done at room temperature however we recommend doing it at 37°C for best consistency with our QC results. Also the TMB incubation time estimate (15-25min) is based on 37°C. |
| Wash buffer | Disolve the wash buffer powder in 1000ml of water to make 1X PBS wash buffer. |
| Biotinylated Anti-Human MIA antibody | It is recommended to prepare this reagent immediately prior to use by diluting the Human MIA Biotinylated antibody (100x) 1:100 with Antibody Diluent. Prepare 100 µl by adding 1 µl of Biotinylated antibody (100x) to 99 µl of Antibody Diluent for each well. Mix gently and thoroughly and use within 2 hours of generation. |
| Avidin-Biotin-Peroxidase Complex | It is recommended to prepare this reagent immediately prior to use by diluting the Avidin-Biotin-Peroxidase Complex (100x) 1:100 with Avidin-Biotin-Peroxidase Diluent. Prepare 100 µl by adding 1 µl of Avidin-Biotin-Peroxidase Complex (100x) to 99 µl of Avidin-Biotin-Peroxidase Diluent for each well. Mix gently and thoroughly and use within 2 hours of generation. |
| Human MIA Standard | It is recommended that the standards be prepared no more than 2 hours prior to performing the experiment. Use one 1ng of lyophilized Human MIA standard for each experiment. Gently spin the vial prior to use. Reconstitute the standard to a stock concentration of 1ng/ml using 1ml of sample diluent. |

| | |
|------------|---|
| | Allow the standard to sit for a minimum of 10 minutes with gentle agitation prior to making dilutions. |
| Microplate | The included microplate is coated with capture antibodies and ready-to-use. It does not require additional washing or blocking. The unused well strips should be sealed and stored in the original packaging. |

Dilution of Human MIA Standard

1. Number tubes 1-8. Final Concentrations to be Tube # 1 – 1000pg/ml, #2 – 500pg/ml, #3 – 250pg/ml, #4 – 125pg/ml, #5 – 62.5pg/ml, #6 – 31.25pg/ml, #7 – 15.625pg/ml, #8 – Sample Diluent serves as the zero standard (0pg/ml).
2. For standard #1, add 1000µl of undiluted standard stock solution to tube #1.
3. Add 300 µl of sample diluent to tubes # 2-7.
4. To generate standard #2, add 300 µl of standard #1 from tube #1 to tube #2 for a final volume of 600 µl. Mix thoroughly.
5. To generate standard #3, add 300 µl of standard #2 from tube #2 to tube #3 for a final volume of 600 µl. Mix thoroughly.
6. Continue the serial dilution for tube #4-7.

Sample Preparation and Storage

These sample collection instructions and storage conditions are intended as a general guideline and the sample stability has not been evaluated.

| Sample Type | Procedure |
|---------------------------|--|
| Cell culture supernatants | Clear sample of particulates by centrifugation, assay immediately or store samples at -20°C. |
| Serum | Use a serum separator tube (SST) and allow serum to clot at room temperature for about four hours. Then, centrifuge for 15 min at approximately 1,000 x g. assay immediately or store samples at -20°C. |
| Plasma | Collect plasma using heparin or EDTA as an anticoagulant. Centrifuge for 15 min at approximately 1,000 x g. Assay immediately or store samples at -20°C. *Note: it is important to not use anticoagulants other than the ones described above to treat plasma for other anticoagulants could block the antibody binding site. |
| Saliva | Collect saliva using a collection device, aliquot and store samples at -20°C. The collection device should not have protein binding or filtering features. |

Sample Dilution

The target protein concentration should be estimated and appropriate sample dilutions should be selected such that the final protein concentration lies near the middle of the linear dynamic range of the assay.

It is recommended to prepare 150 µl of sample for each replicate to be assayed. The samples should be diluted with sample diluent and mixed gently.

Assay protocol

It is recommended that all reagents and materials be equilibrated to 37°C/room temperature prior to the experiment (see Preparation Before The Experiment if you have missed this information).

1. Prepare all reagents and working standards as directed previously.
2. Remove excess microplate strips from the plate frame and seal and store them in the original packaging.
3. Add 100 µl of the standard, samples, or control per well. Add 100 µl of the sample diluent buffer into the zero well. At least two replicates of each standard, sample, or control is recommended.
4. Cover with the plate sealer provided and incubate for 120 minutes at RT (or 90 min. at 37 °C).
5. Remove the cover and discard the liquid in the wells into an appropriate waste receptacle. Invert the plate on the benchtop onto a paper towel and tap the plate to gently blot any remaining liquid. It is recommended that the wells are not allowed to completely dry at any time.
6. Add 100 µl of the prepared 1x Biotinylated Anti-Human MIA antibody to each well.
7. Cover with plate sealer and incubate for 90 minutes at RT (or 60 minutes at 37°C).
8. Wash the plate 3 times with the 1x wash buffer.
 - a. Discard the liquid in the wells into an appropriate waste receptacle. Then, invert the plate on the benchtop onto a paper towel and tap the plate to gently blot any remaining liquid. It is recommended that the wells are not allowed to completely dry at any time.
 - b. Add 300 µl of the 1x wash buffer to each assay well. (For cleaner background incubate for 60 seconds between each wash).
 - c. Repeat steps a-b 2 additional times.
9. Add 100 µl of the prepared 1x Avidin-Biotin-Peroxidase Complex into each well. Cover with the plate sealer provided and incubate for 40 minutes at RT (or 30 minutes at 37°C).
10. Wash the plate 5 times with the 1x wash buffer.
 - a. Discard the liquid in the wells into an appropriate waste receptacle. Then, invert the plate on the benchtop onto a paper towel and tap the plate to gently blot any remaining liquid. It is recommended that the wells are not allowed to completely dry at any time.
 - b. Add 300 µl of the 1x wash buffer to each assay well. (For cleaner background incubate for 60 seconds between each wash).
 - c. Repeat steps a-b 4 additional times.
11. Add 90 µl of Color Developing Reagent to each well. Cover with the plate sealer provided and incubate in the dark for 30 minutes at RT (or 15-25 minutes at 37°C). (The optimal incubation time must be empirically determined. A guideline to look for is blue shading the top four standard wells, while the remaining standards remain clear.)
12. Add 100 µl of Stop Solution to each well. The color should immediately change to yellow.
13. Within 30 minutes of stopping the reaction, the O.D. absorbance should be read with a microplate reader at 450nm.

Data Analysis

Average the duplicate readings for each standard, sample, and control. Subtract the average zero standard O.D. reading.

It is recommended that a standard curve be created using computer software to generate a four parameter logistic (4-PL) curve-fit. A free program capable of generating a four parameter logistic (4-PL) curve-fit can be found online at: www.myassays.com/four-parameter-logistic-curve.assay.

Alternatively, plot the mean absorbance for each standard against the concentration. The measured concentration in the sample can be interpolated by using linear regression of each average relative OD against the standard curve generated using curve fitting software. This will generate an adequate but less precise fit of the data.

For diluted samples, the concentration reading from the standard curve must be multiplied by the dilution factor.

Background on MIA

Melanoma Inhibiting Activity(MIA), also known as cartilage-derived retinoic acid-sensitive protein(CD-RAP), is an approximately 11-15 kDa protein that is expressed as a noncovalent homodimer. MIA is structurally related to OTOX/Otoraplin and MIA-2 in a small family of secreted proteins with one SH3 domain(1-3). And the MIA gene was mapped to 19q13.32-q13.33 by fluorescence in situ hybridization. Besides, It is a marker for malignant melanoma. MIA functions as a chemoattractant for mesenchymal stem cells and enhances their BMP-2 and TGF-beta3 induced differentiation into chondrocytes. MIA-deficient mice exhibit delayed chondrocyte differentiation but enhanced chondrocyte proliferation and cartilage repair.

Submit a product review to Biocompare.com

Submit a review of this product to Biocompare.com to receive a \$20 Amazon giftcard! Your reviews help your fellow scientists make the right decisions. Thank you for your contribution.

