



FREE
VALIDATION



Free Antibody Validation



Experimental Results

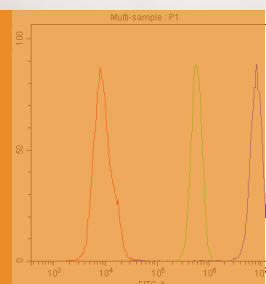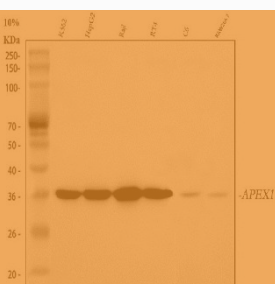

Solve your
validation issues
for free!



Name: Tara MacDonald

Order Number: 2504111036

Lot Number: BS25042391

Antibody Name: Anti-Renalase Rabbit

Monoclonal Antibody

Reviewer: Reviewer: Xiaoqian Cheng


Review Date: 20250429

FREE ANTIBODY VALIDATION

Result Images

IF analysis of Renalase using anti-Renalase antibody (M06291).

Renalase was detected in an immunocytochemical section of mouse Beta-TC6 cells. The cells were fixed with 4% paraformaldehyde for 10 minutes and then treated with a membrane permeabilization agent for 5 minutes. The cells were blocked with 10% goat serum. And then incubated with rabbit anti-Renalase Antibody (M06291) at a dilution of 1:50 overnight at 4°C. DyLight®488 Conjugated Goat Anti-Rabbit IgG (BA1127) was used as secondary antibody at 1:500 dilution and incubated for 30 minutes at 37°C. The section was counterstained with DAPI. Visualize using a fluorescence microscope and filter sets appropriate for the label used.

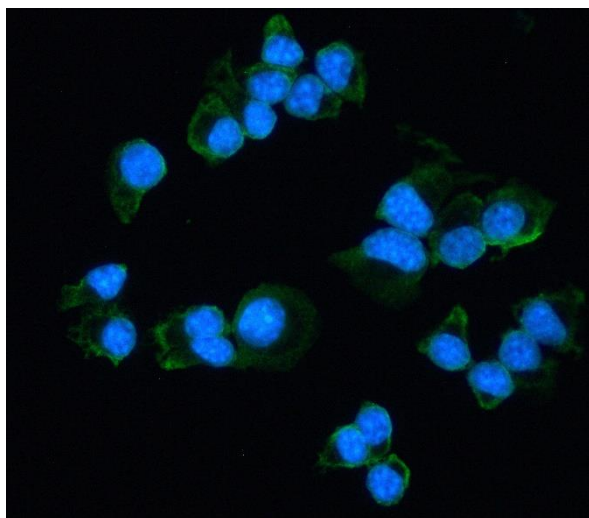

ICC/IF Scoring

★★★★★ (3/5), based on the following criteria:

1. The ICC/IF result showed a secretory expression pattern of Renalase in Beta-TC6 cells, consistent with its predicted secretory characteristics.
2. The ICC/IF result showed a staining pattern consistent with the expected localization.

Evaluation criteria: Experiments rated above ★★ (2/5) are considered a pass, while those rated ★★ (2/5) or less are considered failed.

FREE ANTIBODY VALIDATION

Purchase Recommendation

Corresponding products for the experiment

1. Primary antibody: M06291
2. Secondary antibody: BA1127

If you have any other needs, you may also consider our custom antibody services.

[Have a meeting with us to discuss more experimental details.](#)



FREE ANTIBODY VALIDATION

Summary

Why was this program Created?

Will our antibody, which hasn't been validated on your tissue or cell line, meet your needs? This is a common concern we've heard directly from our customers. In response, we began testing antibodies on unvalidated sample types and providing the results free of charge once they had this question. This practice, born from customer feedback and our dedication to supporting researchers, has consistently received positive recognition.

We formalized it into an official program in 2024 to further honor this commitment. By listening to our customers and going the extra mile, we continue to prioritize your confidence and trust in our products.

Steps (Congratulations! All steps will be free!)

- Enter contact Info
- Find Antibodies
- Select Experiment and samples
- Boster Bio Provides Professional Analysis and Customized Solutions
- Boster Bio Conducts Experiments and Provides You with Reports
- Boster Bio Exclusive Antibody Recommendations
info@bosterbio.com

One-stop solution for your antibody selection needs.



FREE ANTIBODY VALIDATION

Company Promotions



30+years of technique
improvement



20000+ antibodies and
2000+ ELISA kits



60000+ cited
publications



Driven by user's
need



BOSTER
antibody and ELISA experts



(888) 466-3604



www.bosterbio.com



3942 Valley Ave, Pleasanton,
CA 94566, USA



support@bosterbio.com